

SR-5M/6M Stereotaxic Instrument (for Mice)

Stereotaxic Instruments for Mice



SR-5M




SR-6M

These stereotaxic instruments are designed specifically for mice. The mouthpiece, nosepiece, and auxiliary ear bar are all well thought out for improved ease-of-use. The mouthpiece and nosepiece are made slimmer and smaller to fit with the tiny mouth of a mouse. The fore-and-aft and up-and-down adjustments of the mouthpiece/nosepiece have become more sensitive and functional, allowing smooth and quick setup by one-hand maneuvering. The improved auxiliary ear bar allows you to feel the delicate holding touch and the built-in spring permits one-hand maneuvering as well as simple one-touch locking. When not in use, the auxiliary ear bar can be disassembled for cleaning.

Commonly used in conventional models, the 18.70mm AP frame bar permits the use of various accessories other than the supplied SM-15 stereotaxic micromanipulator.

*The easy-to-replace ear bars (pointed tips) of the auxiliary ear bar are available in different angles and diameters. Consult us for specific requirements such as requiring secure fixation or preventing the rupture of an eardrum.

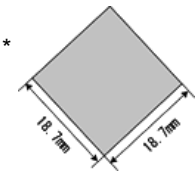

AP frame bar size
Narishige produces world standard 18.7mm square shape AP frame bar.
Please be advised when mounting on your existent frames since our stereotaxic manipulators are tailored for 18.7mm AP frame bar.

Specifications

Accessories	SM-15 Stereotaxic Micromanipulator, EB-5N Auxiliary Mouse Ear Bar
Dimension, Weight	SR-5M: W400 x D300 x H110mm, 8.75kg SR-6M: W400 x D300 x H110mm, 9.2kg